

Invia® Liberty™



SIMPLIFYING NEGATIVE PRESSURE WOUND THERAPY



For User-Friendly Negative Pressure Wound Therapy

Invia® Liberty™


Providing a flexible solution for Negative Pressure Wound Therapy

Negative Pressure Wound Therapy (NPWT) has become a preferred choice of care for treating acute and chronic wounds. Invia® Liberty™ delivers NPWT that satisfies the demands of clinicians and patients.

Invia Liberty's innovative design delivers excellent functionality with its simple set-up and easy operation. Invia Liberty can be used both at hospital and at home, easing transition and facilitating continuity of care.

Patients appreciate the ease-of-use, quiet operation and compact size of the Liberty. The device can be carried in a discreet bag, facilitating mobility. The docking station allows for easy recharging of the Invia Liberty.

Therapy data can be checked on the built-in LCD screen or stored and uploaded to a PC.

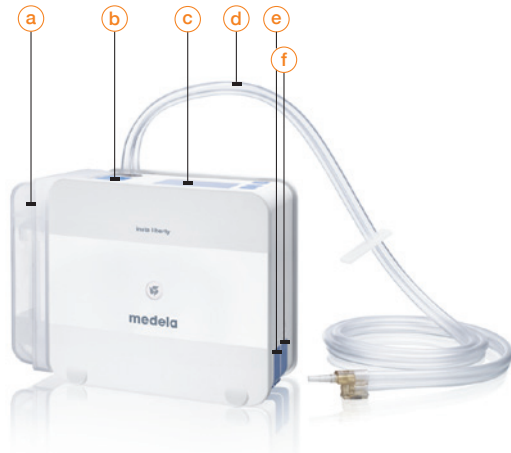

- a Disposable canister 300 ml & 800 ml options with solidifier
- b One-touch canister release button
- c Easy to read LCD screen for intuitive operation and graphic capability
- d Double lumen tubing (PVC) with integrated protection filter and tube clamp
- e Covered, waterspray safe power adapter port
- f Covered, waterspray safe USB port for software updating and therapy history upload to a PC

Indications for Use



- Diabetic/Neuropathic ulcers
- Pressure or venous insufficiency ulcers
- Chronic wounds
- Traumatic wounds
- Acute or sub acute wounds
- Dehisced wounds
- Partial thickness burns
- Flaps and Grafts

Contraindications


- Malignancy in the wound
- Untreated osteomyelitis
- Non-enteric and unexplored fistulas
- Necrotic tissue with eschar present
- Exposed vasculature
- Exposed nerves
- Exposed anastomotic site
- Exposed organs
- Malnourished patients
- Pediatric Use

Medela offers healthcare professionals 50 years of expertise in medical vacuum technology. Our Swiss development and manufacturing process yields mechanical precision to guarantee excellent performance.

Technical data

	- 60 mmHg to 200 mmHg
	Low Flow 5 L/min.

	> 12 Hours
	1 kg 2.2 lbs

Modes:	Constant and Intermittent
Tubing:	Double Lumen Tubing, Y-Connector
Canister:	300 and 800 ml
	H x W x D 95 x 170 x 150 mm 3.74 x 6.69 x 5.91 inch

2 year warranty

Medical Vacuum Technology
for Healthcare Professionals

Please contact us or your local Medela representative for details.

Federal law restricts this device to sale or rental by or on the order of a physician.

Medela and Invia are registered trademarks of Medela Holding AG. Liberty is a trademark of Medela, Inc.

 Medela AG
Lättichstrasse 4b
6341 Baar / Switzerland
www.medela.com

Medela, Inc.
P.O. Box 660
1101 Corporate Drive
McHenry, IL 60050
USA
Phone +1 877 735 1626
Fax +1 815 307 8942
suction@medela.com
www.medelasuction.com

Local contact: