

# Invia® Liberty™



## SIMPLIFYING NEGATIVE PRESSURE WOUND THERAPY



For User-Friendly Negative Pressure Wound Therapy

## Invia® Liberty™


### Providing a flexible solution for Negative Pressure Wound Therapy

Negative Pressure Wound Therapy (NPWT) has become a preferred choice of care for treating acute and chronic wounds. Invia® Liberty™ delivers NPWT that satisfies the demands of clinicians and patients.

Invia Liberty's innovative design delivers excellent functionality with its simple set-up and easy operation. Invia Liberty can be used both at hospital and at home, easing transition and facilitating continuity of care.

Patients appreciate the ease-of-use, quiet operation and compact size of the Liberty. The device can be carried in a discreet bag, facilitating mobility. The docking station allows for easy recharging of the Invia Liberty.

Therapy data can be checked on the built-in LCD screen or stored and uploaded to a PC.

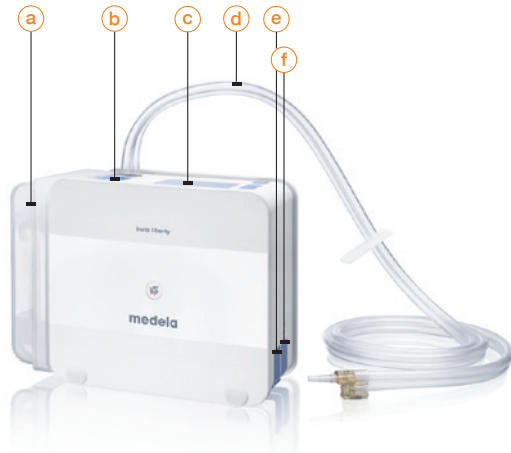

- a Disposable canister 300 ml & 800 ml options with solidifier
- b One-touch canister release button
- c Easy to read LCD screen for intuitive operation and graphic capability
- d Double lumen tubing (PVC) with integrated protection filter and tube clamp
- e Covered, waterspray safe power adapter port
- f Covered, waterspray safe USB port for software updating and therapy history upload to a PC

#### Indications for Use



- Diabetic/Neuropathic ulcers
- Pressure or venous insufficiency ulcers
- Chronic wounds
- Traumatic wounds
- Acute or sub acute wounds
- Dehisced wounds
- Partial thickness burns
- Flaps and Grafts

#### Contraindications


- Malignancy in the wound
- Untreated osteomyelitis
- Non-enteric and unexplored fistulas
- Necrotic tissue with eschar present
- Exposed vasculature
- Exposed nerves
- Exposed anastomotic site
- Exposed organs
- Malnourished patients
- Pediatric Use

Medela offers healthcare professionals 50 years of expertise in medical vacuum technology. Our Swiss development and manufacturing process yields mechanical precision to guarantee excellent performance.

#### Technical data

	- 60 mmHg to 200 mmHg
	Low Flow 5 L/min.

	> 12 Hours
	1 kg 2.2 lbs

Modes:	Constant and Intermittent
Tubing:	Double Lumen Tubing, Y-Connector
Canister:	300 and 800 ml
	H x W x D 95 x 170 x 150 mm 3.74 x 6.69 x 5.91 inch

2 year warranty

Medical Vacuum Technology  
for Healthcare Professionals

**Please contact us or your local Medela representative for details.**

Federal law restricts this device to sale or rental by or on the order of a physician.

Medela and Invia are registered trademarks of Medela Holding AG. Liberty is a trademark of Medela, Inc.

 Medela AG  
Lättichstrasse 4b  
6341 Baar / Switzerland  
www.medela.com

Medela, Inc.  
P.O. Box 660  
1101 Corporate Drive  
McHenry, IL 60050  
USA  
Phone +1 877 735 1626  
Fax +1 815 307 8942  
suction@medela.com  
www.medelasuction.com

#### Local contact: